

# Miru

## 1month Menicon

Silicone hydrogel monthly  
disposable lenses

A totally breathable lens that ensures  
every part of the eye gets the oxygen it  
needs for healthy, happy eyes.



# Technology in balance



Health



Vision



Comfort

Miru 1 month, a unique technology designed to meet the demands of today's contact lens wearer.

---

## Material and surface technologies

---

### MeniSilk™

Unique silicone polymerization, innovative hydrophilic monomer

- Ultra high Dk/t
- Exceptional hydration
- Optimized transparency


### NanoGloss™

Unique surface technology, nanometer precision

- Super smooth surface
- Resistance to bacteria
- Excellent wettability

MeniSilk™

NanoGloss™

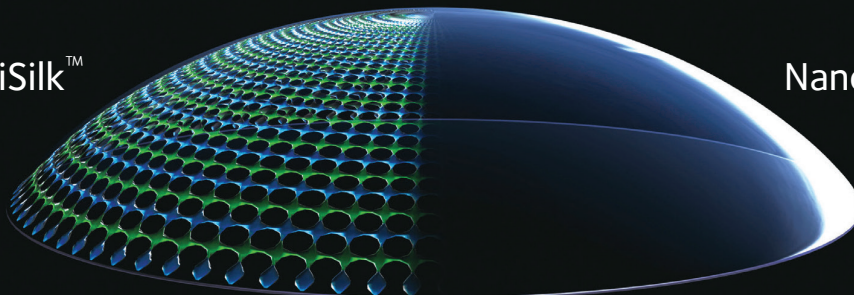


---

# Miru 1 month toric


---

## Visiostable design™

Horizontal dynamic  
stabilization zones

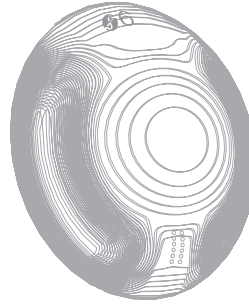


Unique double vertical  
asymmetric Thin Zones

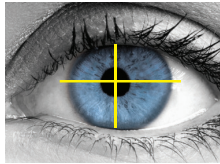

## Anatomical profile


The unique asymmetric vertical Thin Zones match the eyelids' natural asymmetric coverage of the cornea, harnessing the natural lid force, optimizing centration and preventing rotation.

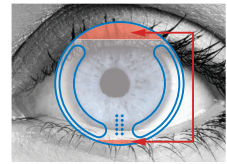

---

# Miru 1 month multifocal

---

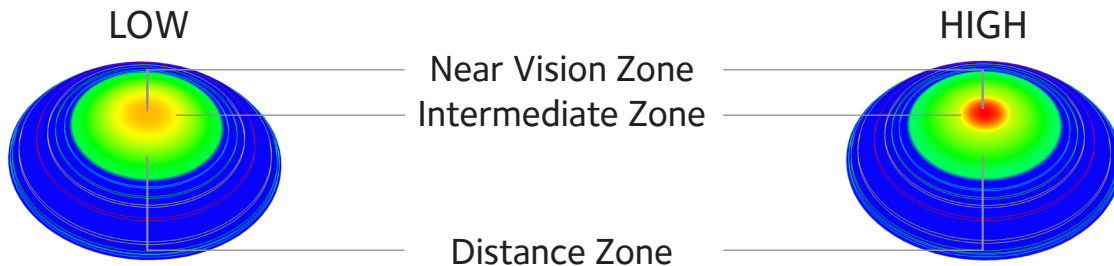

## Wearer Benefit

The 3 zones are designed in optimal balance to control depth of focus, providing a natural transition between near and distance vision.

<div>Miru</div> <div>1month Menicon</div>		Product Specifications		
		Sphere	Toric	Multifocal
Characteristics	Material	asmoofilcon A (Silicone hydrogel)		
	Water Content	40%		
	Dk/t @ -3.00D	161x 10 <sup>-9</sup> (cm/sec) • (mLO <sub>2</sub> /(mL x mmHg))		
	Center Thickness	0.08mm @ -3.00D		
Parameters	Diameter	14.00mm		14.20mm
	Base Curve	8.3mm/8.6mm	8.6mm	8.6mm
	Sphere	+6.00D to -13.00D	plano to -10.00D	+6.00D to -13.00D
		<div> <div>*No plano</div> <div> <div>0.25D Steps</div> <div>0.50D Steps</div> </div> <div> <div>+6.00</div> <div>-6.00</div> <div>-13.00</div> </div> </div>		<div> <div>0.25D Steps</div> <div>0.50D Steps</div> </div> <div> <div>+6.00</div> <div>-6.00</div> <div>-13.00</div> </div>
		<div> <div>0.25D Steps</div> <div>0.50D Steps</div> </div> <div> <div>plano</div> <div>-6.00</div> <div>-10.00</div> </div>		<div> <div>0.25D Steps</div> <div>0.50D Steps</div> </div> <div> <div>+6.00</div> <div>-6.00</div> <div>-13.00</div> </div>
	Cylinder		-0.75, -1.25, -1.75	
	Axis		10° to 180° (in 10° steps)	
	Addition			LOW(+1.00), HIGH(+2.00)
Lens Marking				
Wear		Daily wear monthly replacement		
Packaging		Available in 6 packs		